


microtube

- IDEAL FOR LOW TEMPERATURE STORAGE
- SUITABLE FOR 0.5ML AND 1.5ML MICROTUBES + 2.0ML CRYOTUBES
- IDEAL LOW COST TRANSPORTATION
- COMPACT DESIGN
– 185 x 110 x 60MM

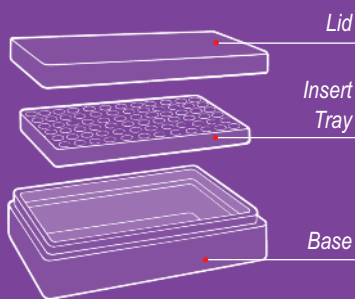

EPRAK STORAGE

The Eprak is a low cost multi-purpose box manufactured in expanded polystyrene.

With its removable Insert Tray, it is ideal for use as both a work station and storage box for microcentrifuge tubes, of both 1.5ml (50) and 0.5ml (72) capacities as well as 2.0ml cryotubes (36).

Eprak can be used along with either liquid nitrogen or dry ice for flash freezing of samples, or the expanded polystyrene Insert Tray can be floated in a thermostatically controlled bath to allow defrosting or incubation of tube contents.




- Expanded polystyrene allows controlled temperature changes
- Can be used with crushed or dry ice, water and liquid nitrogen
- Removable Insert Tray floats in water bath
- Defrosts or incubates tube contents
- Inserts can easily be trimmed to a smaller size

coolcube

- MAINTAINS 0°C FOR 4 HOURS
- HOLDS BOTH MICROTUBES AND A PCR PLATE
- CAN BE USED AS A PCR WORKSTATION

CoolCube is designed to store samples safely on the lab bench without sample degradation caused by temperature increase or fluctuation. CoolCube maintains temperature at approximately 0°C for up to four hours. Simply place it in a freezer overnight and it will be ready to use.

One side of the CoolCube can be used as a convenient PCR workstation for 12 standard 1.5/2.0ml microtubes with one 96 well PCR plate, or with 12x 0.2ml strips or 96x0.2ml tubes. The other side may be used to hold up to 36 standard 1.5/2.0ml microtubes.

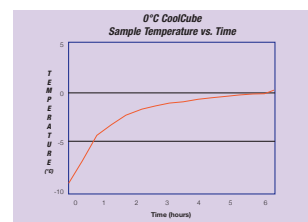

Technical Specifications

TEMPERATURE:	0°C FOR UP TO 4 HOURS
CAPACITY:	36 x 1.5ML 96 x 0.2ML 12 x 0.2ML PCR STRIPS 1 x 96 WELL PCR PLATE
DIMENSIONS:	142 x 142 x 63.5MM
WEIGHT:	900G
MATERIAL:	DURASTAR POLYMER

ICECOOL baths



- DURABLE PLASTIC
- EASILY CLEANED
- EXCELLENT LONG TERM THERMAL STORAGE

IceCool Baths are constructed of durable plastic and offer both robustness and long life. They will also withstand reasonable use with dry ice. Lids are available for the large bucket.

ORDERING INFORMATION

BC-5003	EPRAK with 1.5ml Microtube Inserts, pk/30
BC-5004	EPRAK with 0.5ml Microtube Inserts, pk/30
BC-5008	EPRAK with 2.0ml Cryotube Inserts, pk/30

BC-5005	Replacement 1.5ml Insert Trays, pk/30
BC-5006	Replacement 0.5ml Insert Trays, pk/30
CLCB	Coolcube Storage Rack

CSLIB	Large IceCool Bath - White	CSLIBLID	Lid for Large IceCool Bath - White	CSMIB	Mini IceCool Bath - White
-------	----------------------------	----------	------------------------------------	-------	---------------------------