

microcentrifuges

- COMBI-ROTOR FOR TUBES AND STRIPS
- DOUBLE THE CAPACITY OF TRADITIONAL MINI CENTRIFUGES
- NEARLY SILENT OPERATION
- STARTS & STOPS WITH OPENING/CLOSING OF THE LID

MULTIFUGE




Unlike traditional mini centrifuges, the new multiFuge 12 eliminates the need to change rotors when switching between microtubes and PCR strips. The included, unique COMBI-Rotor is all that is required for running 12 microtubes and 4 PCR strips simultaneously. With a fixed speed that produces 2,000 x g, this centrifuge is perfect for quick spin downs. Simply close the lid and the unit quickly ramps up to 5500 rpm. Open the lid, and the rotor quickly decelerates for removal of samples.

At just 14cm wide and less than 11cm high, the multiFuge 12 truly is a personal centrifuge with unmatched capacity and flexibility.

Technical Specifications

SPEED:	5,500 RPM / 2,000 x G
CAPACITY:	12 x 1.5 / 2.0 ML TUBES, 32 x 0.2 ML PCR TUBES, 4 x PCR STRIPS (8x0.2 ML)
DIMENSIONS:	14 x 20 x 11.2 CM
WEIGHT:	5 KG

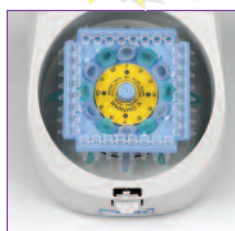

Unique COMBI-Rotor



Compact, low profile design

QuickSpin

- SUPPLIED WITH BOTH STANDARD MICROTUBE AND STRIP TUBE ROTORS
- IDEAL FOR QUICK SPIN-DOWNS AND MICROFILTRATION
- STARTS AND STOPS IN SECONDS
- COMPACT DESIGN

Quickspin is perfect for microfiltration and rapid spin-down of sample from the walls and caps of microcentrifuge tubes.

Occupying less than 6 inches square of bench space, the Quickspin has a very small footprint, making it easy to use in the lab. Rotors and adaptors, which may be easily interchanged, are supplied as standard to accommodate 1.5 ml, 0.5 ml and 0.4 ml tubes, as well as 0.2 ml strips and individual tubes.

A highly durable stainless steel hinge pin facilitates easy opening of the translucent lid, while an on/off switch is located on the side of the centrifuge to start and stop operation. Alternatively, with the switch in the 'on' position, the centrifuge can be started and stopped by closing and opening the lid.



Technical Specifications

MAXIMUM SPEED	6000RPM
MAXIMUM G FORCE	2,000 x G
CAPACITY	6 x 1.5/2.0 ML 2 x 0.2ML STRIPS
DIMENSIONS (W x D x H)	15 x 15 x 11.7 CM
WEIGHT	0.45 KG

ORDERING INFORMATION

MFUGE	MultiFUGE™ Mini Centrifuge, with COMBI-Rotor for 12 microtubes and PCR strips/tubes	C0812-A5	Adapters, 0.5 ml, pack of 6
		C0812-A2	Adapters, 0.2 ml, pack of 6
QSPIN	QuickSpin complete with 1.5/2.0 ml rotor, strip tube rotor, 0.5 and 0.4 ml adapters, Purple lid		